

# CRImage

## a package for classifying cells and calculating tumour cellularity

Henrik Failmezger\*, Yinyin Yuan, Florian Markowetz

\* E-mail: failmezger@cip.ifi.lmu.de

## Contents

1	Load the package	1
2	Threshold an image	1
3	Segment an image	2
4	Creating a training set	2
5	Creating the classifier	4
6	Classification of cells	4
7	Calculation of cellularity	5
8	Classification of Aperio Image Slides	5

## 1 Load the package

The package is loaded by the following command:


```
> library(CRImage)
```

See the functions of the package EBIImage to read, write and manipulate an image.

## 2 Threshold an image

To create a binary image of a coloured image the function `calculateThreshold` can be used. This function calculates Otsu threshold of a grayscale image.

```
> f = system.file("extdata", "exImg2.jpg", package = "CRImage")
> exImgB = readImage(f)
> indexWhitePixel = which(exImgB[, , 1] > 0.85 & exImgB[, , 2] >
+ 0.85 & exImgB[, , 3] > 0.85)
> exImgB = channel(exImgB, "gray")
> t = calculateThreshold(as.vector(exImgB[-indexWhitePixel]))
> exImgB[which(exImgB >= t)] = 1
> exImgB[which(exImgB < t)] = 0
```

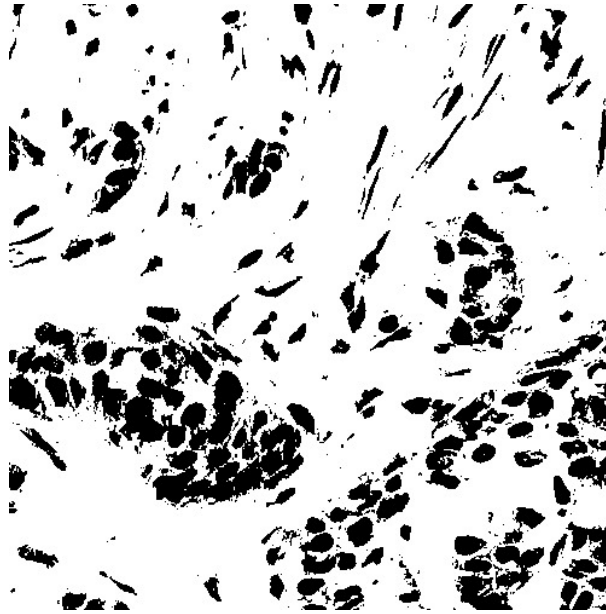

**Figure 1.** Thresholded image. Cell nuclei are coloured white, background is coloured black.

The image is read with `readImage`. `indexWhitePixel` finds white pixels in the image. These pixels are considered as background and excluded from thresholding. The image is converted to grayscale by the function `channel(exImgB,"gray")` of `EImage` and thresholded with `calculateThreshold`. The binary image is created by finding pixels with gray values larger or smaller than the threshold and colouring them black or white.

### 3 Segment an image


An image can be segmented to find cells in the image, using the function `segmentImage`.

```
> f = system.file("extdata", "exImg2.jpg", package = "CRImage")
> segmentationValues = segmentImage(f, maxShape = 800, minShape = 40,
+   failureRegion = 2000)
```

The image is converted to grayscale and thresholded. Morphological opening is used to delete clutter and to smooth the shapes of the cells. The watershed algorithm is used to separate clustered cells. The parameter `maxShape` defines the maximal shape of cell nuclei. Segmented nuclei which exceed this value will be thresholded and segmented again. The parameter `minShape` defines the minimum size of cell nuclei. Cell nuclei, which fall below this value will be deleted. The parameter `failureRegion` defines, when artifacts in the image should be deleted. Dark regions which exceed this value will be deleted. The list element `segmentationValues[[1]]` holds the original image, `segmentationValues[[2]]` holds the segmented image and `segmentationValues[[3]]` holds features, which were calculated for the segmented objects. The segmented objects can be drawn by the function `display` (see `EImage`).

### 4 Creating a training set

To classify cells in images first an appropriate training set has to be created.

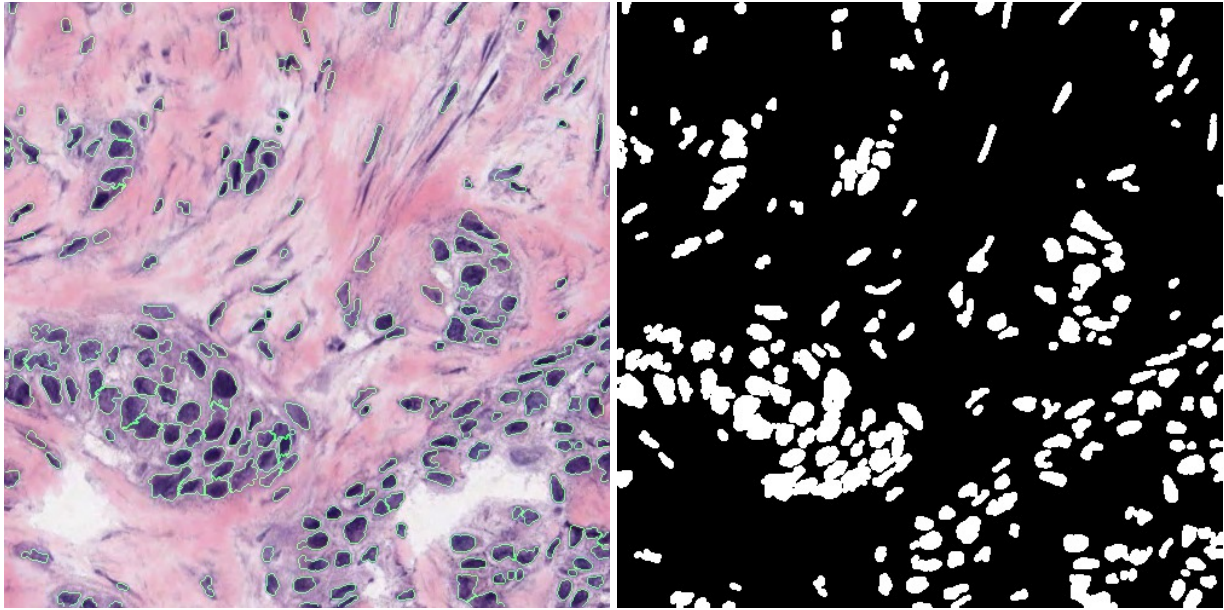

**Figure 2.** SegmentedImage. Cell nuclei are labeled in green.

```
f = system.file("extdata", "exImg.jpg", package="CRImage")
trainingValues=createTrainingSet(filename=f,maxShape=800,minShape=40,failureRegion=2000)
```

The list trainingValues returns two values. The first value is an image in which every segmented cell is numbered. The second value is a table with features for every cell:


index	class	g.x	g.y	g.s	g.p	g.pdm	g.pdsd	g.effr...
1	<NA>	148.2203	4.855932	118	36	5.721049	1.0130550	6.128668...
2	<NA>	160.2719	4.763158	114	38	5.659780	1.0331399	6.023896...
3	<NA>	183.7975	3.101266	79	35	4.593585	0.9575905	5.014627...
4	<NA>	196.3500	4.242857	140	43	6.424591	1.4433233	6.675581...
5	<NA>	271.5694	2.680556	72	29	4.504704	1.2490794	4.787307...
6	<NA>	338.5221	6.530973	113	35	5.601531	1.0849651	5.997418...
7	<NA>	456.0179	2.946429	112	39	5.726556	1.8101368	5.970821...
8	<NA>	556.3778	9.018519	270	81	9.196894	2.4687870	9.270581...
9	<NA>	575.1777	10.935950	484	101	11.959995	2.1378649	12.412171...
10	<NA>	592.7391	1.724638	69	41	5.145094	2.2782539	4.686511...

To create the training set, class values for the cells have to be inserted in the column "class":

index	class	g.x	g.y	g.s	g.p	g.pdm	g.pdsd	g.effr
1	normal	148.2203	4.855932	118	36	5.721049	1.0130550	6.128668
2	malignant	160.2719	4.763158	114	38	5.659780	1.0331399	6.023896

The values in the column "index" match the numbers for the cells in the image. You can save the table as tab separated by:

```
write.table(trainingData,file="path",sep="\t",rownames=F)
```

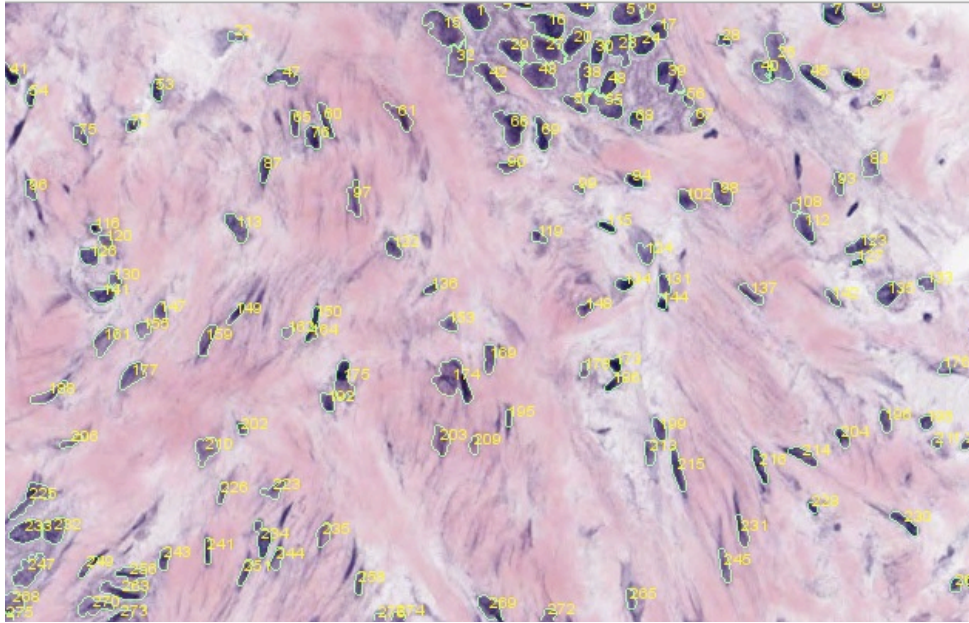

**Figure 3.** The image with labeled cell nuclei. The number of the cell nuclei correspond to the index in the feature table.

and open it for example in a spreadsheet program, however be careful not to shift the columns. You do not have to assign a class value to every cell, but do ensure that there are enough examples for every class. Class values can be numbers (e.g. 1,2,3) or strings (e.g. "normal", "malignant"). You can choose at most 10 classes. If you want to use the kernel smoothing approach for classification you can only choose two classes or you have to specify a cancer class, when classifying cells.

## 5 Creating the classifier

The command:


```
> f = system.file("extdata", "trainingData.txt", package = "CRImage")
> trainingData = read.table(f, header = TRUE)
> classifierValues = createClassifier(trainingData, topo = FALSE)
> classifier = classifierValues[[1]]
```

creates the classifier. The classifier is a Support Vector Machine provided by the package e1017.

## 6 Classification of cells

After having created the classifier, cells in another image can be classified (Figure 2).

```
> f = system.file("extdata", "exImg2.jpg", package = "CRImage")
> classValues = classifyCells(classifier, filename = f, KS = TRUE,
+   maxShape = 800, minShape = 40, failureRegion = 2000)
```

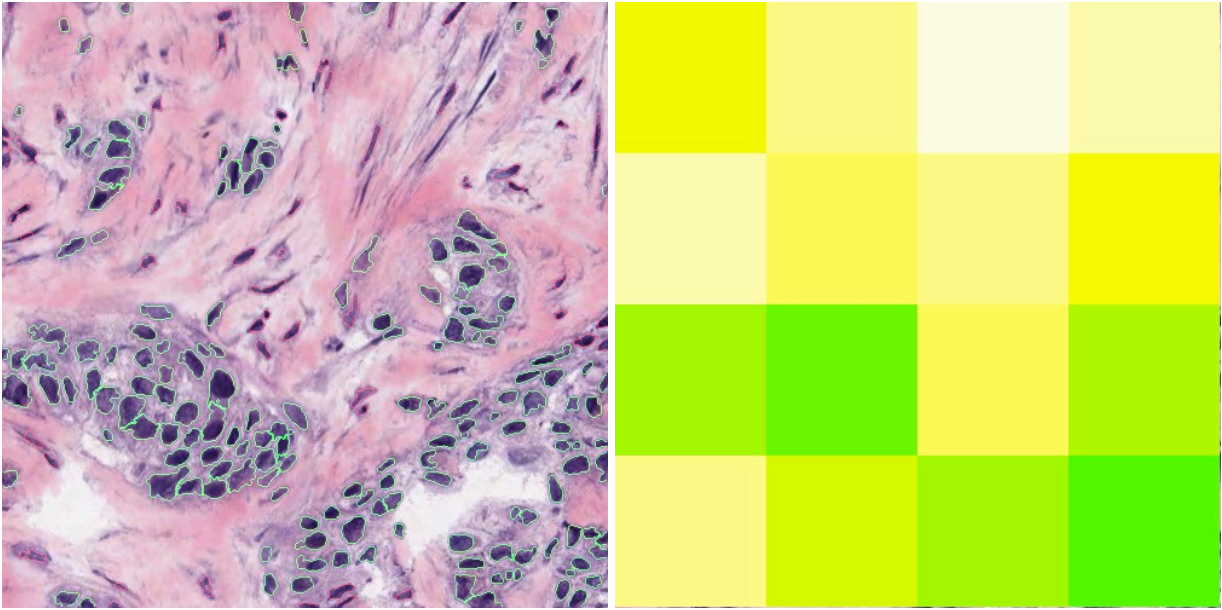

**Figure 4.** Classified Image (left) and heatmap of cellularity values (right). Malignant cells are colored green, whereas other cells are coloured red. The heatmap of tumour cellularity indicates regions of large tumour cellularity (strong green values) and low tumour cellularity (white regions).

## 7 Calculation of cellularity

If tumour images are processed, the cellularity of the tumour can be calculated. The image is first segmented and the cell types are classified. Afterwards, the cellularity of the tumour is determined. The function needs to know, which class value should be the tumour class.

```
> t = system.file("extdata", "trainingData.txt", package = "CRImage")
> trainingData = read.table(t, header = TRUE)
> classifier = createClassifier(trainingData, topo = FALSE)[[1]]
> f = system.file("extdata", "exImg2.jpg", package = "CRImage")
> exImg = readImage(f)
> cellularityValues = calculateCellularity(f, classifier = classifier,
+   cancerIdentifier = "1")
```

## 8 Classification of Aperio Image Slides

Large pathological images are often difficult to process due to their often large file size. Images of ScanScope TX scanner can be saved in the CWS file format. In this format the images are separated in smaller subimages. These images can be process with the function `processAperio`.

```
> t = system.file("extdata", "trainingData.txt", package = "CRImage")
> trainingData = read.table(t, header = TRUE)
> classifier = createClassifier(trainingData, topo = FALSE)[[1]]

> dir.create("AperiOutput")
> f = system.file("extdata", package = "CRImage")
```

```

> f = file.path(f, "8905")
> processAperio(classifier = classifier, inputFolder = f, outputFolder = "AperiOutput",
+   identifier = "Da", numSections = 2, cancerIdentifier = "c",
+   maxShape = 800, minShape = 40, failureRegion = 2000)

```

The function gets a classifier, which was for instance created with `createClassifier`. An input folder has to be specified, which includes the subimages. The common identifier of the subimages has to be specified (either "Da" or "Ss", depending on the adjustments of the Aperio software). If the image includes different sections for which cellularity has to be calculated separately, this can be specified with `numSections`, the function will afterwards do an automatical clustering of the subimages to their corresponding sections. The function needs to know, which class of the classifier is the tumour class, which is specified with `cancerIdentifier`. The function creates three folders in the output folder. The folder `classifiedImage` includes the classified subimages. The folder `Files` includes files with cellularity values for every section. The folder `tumourDensity` includes cancer heatmaps for every subimage.